

96-Well PCR Plates Non-shirted

For medium and high sample throughput, 96-well plates from Medicalet are the efficient solution for PCR and qPCR and can be conveniently filled with multichannel pipettes or pipetting robots.

The uniform and thin wall thicknesses of the wells ensure rapid transfer of the temperature from the cycler to the sample, thus reducing cycle times.

A number of skirt choices ensures a good fit in common thermal cyclers for efficient heat transfer. For high flexibility, unskirted plates are the first choice.

- ✓ Broad thermal cycler compatibility.
- ✓ High-contrast, black alphanumeric indexing for easy well identification.
- ✓ Crystal-clear wells to easily verify sample volume.
- ✓ Ultra-thin, uniform wells ensure optimal heat transfer for high reaction efficiency.
- ✓ Well ridges reduce sample evaporation when sealing with film, foil, or strip caps.
- ✓ Bar-coded semi-skirted and full-skirted plates available.
- ✓ Manufactured from high-clarity, ultra-pure polypropylene (meets FDA 21 CFR 177.1520 & USP Class VI).
- ✓ Certified free of detectable DNase, RNase, DNA, PCR inhibitors, and tested pyrogen-free.

Applications

- ✓ PCR
- ✓ qPCR/Real-Time PCR
- ✓ Sequencing


Order Information

96-Well PCR Plates Non-shirted, free of DNase, RNase & DNA

Cat. No.	Product	Capacity	color	Type	Packing	Qty. per box/case
PCR022096	96-Well PCR Plates	0.2 mL	White	Non-shirted	Box	5 / 10
PCR022196	96-Well PCR Plates	0.2 mL	Colorless	Non-shirted	Box	5 / 10
PCR012096	96-Well PCR Plates	0.1 mL	White	Non-shirted	Box	5 / 10
PCR012196	96-Well PCR Plates	0.1 mL	Colorless	Non-shirted	Box	5 / 10



Technical Specification

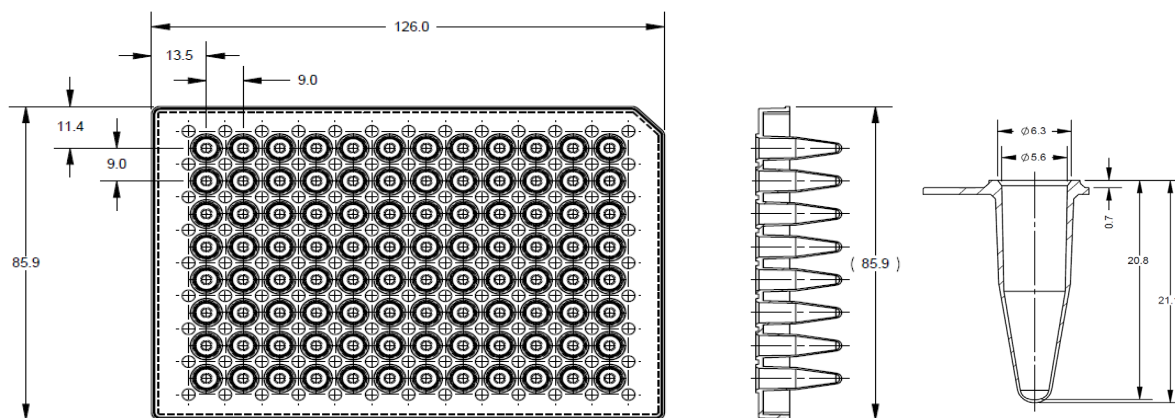

96-Well PCR Plates Non-shirted

96-Well PCR Plates Semi-shirted


Order Information

96-Well PCR Plates Semi-shirted, free of DNase, RNase & DNA



Cat. No.	Product	Capacity	color	Type	Packing	Qty. per box/case
PCR021096	96-Well PCR Plates	0.2 mL	White	Semi-shirted	Box	5 / 10
PCR021196	96-Well PCR Plates	0.2 mL	Colorless	Semi-shirted	Box	5 / 10
PCR011096	96-Well PCR Plates	0.1 mL	White	Semi-shirted	Box	5 / 10
PCR011196	96-Well PCR Plates	0.1 mL	Colorless	Semi-shirted	Box	5 / 10

Technical Specification

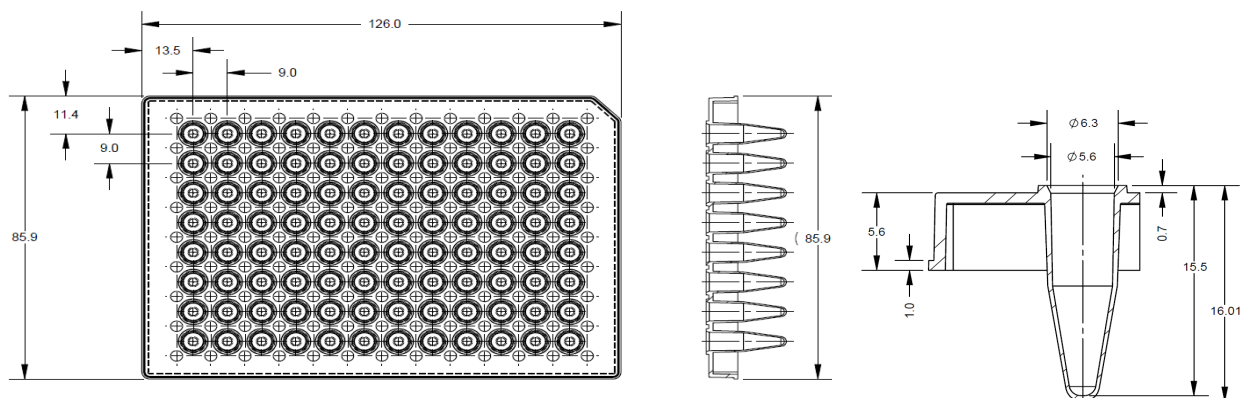

96-Well PCR Plates Semi-shirted